


AaRiss™

Drug Eluting Balloon

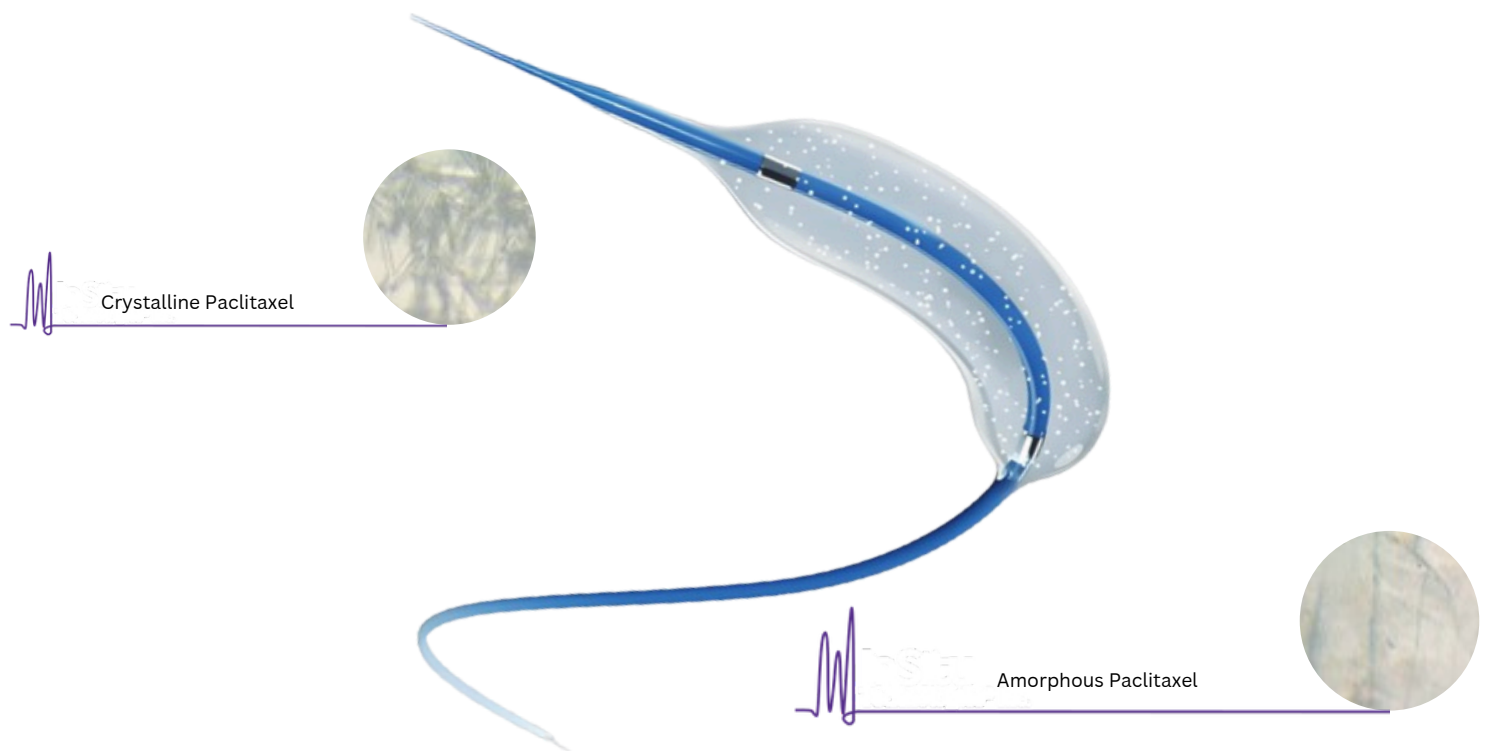


- Clinical evidence shows a low rate of restenosis
- The unique paclitaxel coating technology reduces drug
- Non-expedition design
- Different drug dose can be used for coronary artery stenosis ($3\mu\text{g}/\text{mm}^2$) and intracranial artery stenosis ($2\mu\text{g}/\text{mm}^2$)

COATING TECHNOLOGY

Optimized coating formulation

The PurePac™ coating technology provides an optimal delivery mechanism, ensuring minimal drug loss during tracking and inflation, rapid drug transfer to the vessel wall and a sustained therapeutic effect.

Drug	Paclitaxel anti-proliferative agent
Dose	Coronary artery stenosis (3 µg/mm²) and intracranial artery stenosis (2 µg/mm²)
Coating	Pure paclitaxel, non-ecipient design
Layers	Composite coating, 70% crystalline and 30% amorphous

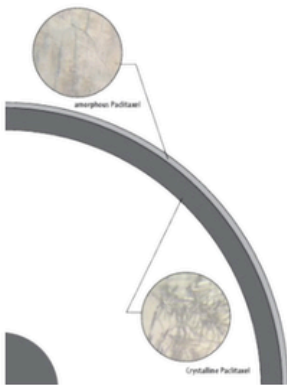

Technical Data

Specifications	Description
Balloon System	Rapid Exchange
Balloon Material	Semi-compliant Polymer (Paclitaxel Coated)
Balloon Fold	Tri-Fold Design
Compliance	Semi-compliant
Markerbands	Two marker bands
**Nominal Pressure	6 atm
Rated Burst Pressure	14 atm
Shaft Size (Proximal/Distal)	Proximal: 0.68 ± 0.05 mm Distal: 0.90 ± 0.10 mm
*Minimum Guiding Catheter	5F
Tip Entry Profile	Low profile soft tapered tip for smooth lesion entry
Recommended Guidewire	0.014"
Delivery Catheter Length	142 ± 4 cm
Depth Markings	for accurate positioning under fluoroscopy
Balloon Diameter	2.00-4.00 mm
Balloon Length	6mm-40mm

* See individual manufacturer for Guide Catheter Fr equivalent
** In Vitro testing data, Balloon Catheter may be affected by lesion characteristics.

AaRiss™ Ordering Information

Balloon Diameter (mm)	Balloon Length (mm)			
2.0	06, 08, 10, 12, 15, 20, 35, 40	92	Balloon Diameter	Balloon Length
2.25	06, 08, 10, 12, 15, 20, 25, 30, 35, 40			
2.5	06, 08, 10, 12, 15, 20, 25, 30, 35, 40			
2.75	06, 08, 10, 12, 15, 20, 25, 30, 35, 40			
3.0	06, 08, 10, 12, 15, 20, 25, 30, 35, 40			
3.25	06, 08, 10, 12, 15, 20, 25, 30, 35, 40			
3.5	06, 08, 10, 12, 15, 20, 25, 30, 35, 40			
3.75	06, 08, 10, 12, 15, 20, 25, 30, 35, 40			
4.0	06, 08, 10, 12, 15, 20, 25, 30, 35, 40			

The content in this brochure is provided for informational purpose only and is not intended or recommended as a substitute for educational professional medical advice or diagnosis purposes. The information is for use only in countries with applicable health authority product registration.

FDA Registered Facility
EN ISO 13485:2016 Certified
For International Distribution. Rev. 11/15/2025
©2025 COPYRIGHT InSitu Technologies® Inc.

Contact Us

1301 Corporate Center Drive, Suite 160,
Eagan, 55121, Minnesota, United States of
America

hello@insitu-tech.com

651-389-1017